

# **Sung Kwang Meditech**

# For surgical patient & Rehabilitation Radiology Products

# **Operating Room**

Snow-pac Several device Position pad Gel product Case Silicon product



# Rehabilitation


Gel wheelchair cushion pad Gel wheelchair back cushion pad



Silver
Gel embo cushion mattress





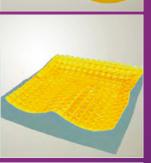

# **Radiation Protection For:**

CardiacCathLab InterventionalRadiology Pain Management











# **Operating Room**

# SNOW PAC


# FOR SURGERY POSITIONING PATIENTS

**Snow Pac** provides more stable and comfortable positioning for patient, saves valuable time for the entire O.R. team in every case.

It prevents pressure points. It aids faster, easier set-up.

aids to achieve comfortable position. aids better surgery position presentation.

# PRINCIPLE OF SYSTEM

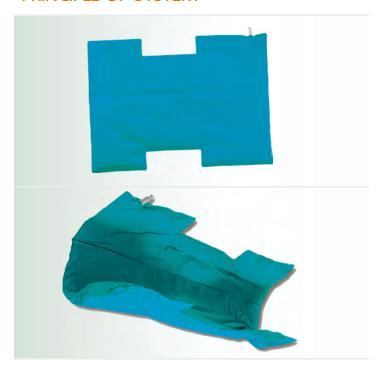

#### **Before**

Snow-Pacs are soft pads (available in various sizes and shapes) that are filled with thousands of tiny plastic beads.

The Snow-Pac is molded around the patient and then a suction line is attached to it.

As a vacuum is created inside the Snow-Pac within 15 seconds, outside atmospheric pressure forces the beads together so they cannot move.

#### **After**

When the suction line is disconnected, the Snow-Pac retains its firm, molded shape and holds the patient securely in position. After surgery, the valve is opened, and as air reenters into the Snow-Pac it becomes soft again, ready for reuse.

### **FFFFCT**

Because Snow Pacs are molded completely and firmly to shape of the patient's body and more securely than sandbags, towel rolls and other devices that are merely propped against the patient.

#### More Comfortable

The apparent wrinkles are inverted-drawn inward by suction and will not irritate the skin.

#### BetterPrese ntation

With Snow Pacs the patient can be positioned exactly as the surgeon desires for the best presentation of the operative site. Snow Pacs conform to any shape, so the patient's position does not have to be compromised.

#### More Versatile

With Snow Pacs patient can often be placed in positions that could not be easily achieved using other positioning devices

# NoPressu re Points

A Snow-Pac supports all parts of the patient's body uniformly, like a floatation mattress Weight is distributed evenly, instead of being supported only by the bony prominences, which Leads to pressure point injuries.

Patients can be positioned much more quickly and easily using Snow Pacs, and if necessary, be repositioned during a case.

### OTHER FEATURES

### Any Vacuum Source

Use wall suction, aspirator pump or any other convenient O.R. vacuum source.

### Reusable

Snow Pacs give years of O.R. service. Rubberized vinyl fabric is soft but durable.

### Radiolucent

Excellent radio transparency with only minimal and clearly identifiable artifacts.

#### Clean, Disinf ect, Sterilize

Can be cleaned with regular soaps and liquid disinfectants or sterilized by gas.

### Saves Money

Snow Pacs save the cost of sandbags, towel rolls and other positioning devices, and they save time in every case.

#### Repair Service

Factory inspection and repair service is available at nominal cost after warranty period Full One Year Warranty.

# Sung Kwang Meditech < €

# **INDICATION**

# KIDNEY POSITION, LATERAL POSITION







A-1 type A-1-1 type



A-1-2 type A-1-3 type

# Easily up in nephrectomy (A type), and also for handicapped person with cerebral palsy.(A-1 type).



A type	73(L)* 100(W)cm
A-1 type	73(L) * 100(W)cm
A-1-1 type	90(L) * 100(W)cm
A-1-2 type	73(L) * 100(W)cm(U)
A-1-3 type	90(L) * 100(W)cm(U)

Provides hyperextension, prevents postoperative discomfort.

# **THYROID POSITION**




B type B-1 type

**B** type 72(L) \* 88(W)cm B-1 type 118(L) \* 88(W)cm

Stabilizes head in eye, facial or dental surgery.

# **HEAD STABILIZATION POSITION**



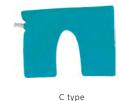

**SUPINE POSITION** 



**VARIOUS FULL BODY POSITION** 








D type

C type 30(L) \* 40(W)cm

Supports and eliminates back pains during long procedures. Prevents foot drop splint in Emergency.



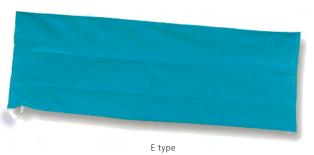

E type 165(L) \* 72(W)cm

# **Operating room**

# **CASES**




inner: magnet For surgery in OR It aids to count needle thoroughly.



SK-18

### **NEEDLE HUNTER**

It finds needles by rolling on the floor

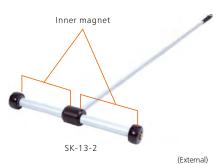

17(L) \*8 (w) \* 2(H)cm 14(L) \* 4(w) \* 1.5(H)cm 100(L) \* 30(w) SK-18-1 (L) SK-18 (S) SK-13-2

### MICRO INSTRUMENT ONE FLOOR CASE





SK-19-a(1F)

26(L) \* 17(w) \* 2.5(H)cm

**MICRO INSTRUMENT** TWO FLOOR CASE



SK-19 (2F)



SK-19 (2F)

For microsurgical and other delicate instruments, for maximum protection and flexible handling, made of resilient polymer. Protruded silicon embomat is placed on the bottom prevents instrument to move freely. Cases are designed to withstand repeated exposure to 134°C steam temperature.

Steam / EO / Plasma Sterilizable

(External)

SK-19(2F)

26(L) \* 17(W) \* 4(H)cm

# Sung Kwang Meditech ← €

# MICRO INSTRUMENT ONE FLOOR CASE





For microsurgical and other delicate instruments, for maximum protection and flexible handling, made of resilient polymer.

Protruded silicon embomat is placed on the bottom prevents instrument to move freely.

Cases are designed to withstand repeated exposure to 134°C steam temperature.



Steam / EO / Plasma Sterilizable




(Extornal)

SK-19-s	16(L) * 6(w) * 3(H) cm
SK-19-1	25(L) * 5 (W) * 2.5(H)cm
SK-19-3	30(L) * 10(W) * 2.5(H)cm

# **Operating room**

### MICRO INSTRUMENT ONE FLOOR CASE



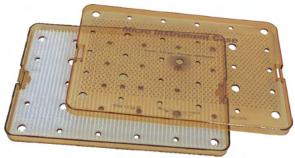

SK-19-5-a



# Steam / EO / Plasma Sterilizable

(External)

SK-19-5-a(1F) 38(L) \* 26(w) \* 2.5(H)cm SK-19-5-b(1F) 39(L) \* 27(w) \* 12(H) cm

# MICRO INSTRUMENT TWO FLOOR CASE



For microsurgical and other delicate instruments, for maximum protection and flexible handling, made of resilient polymer.
Protruded silicon embomat is placed on the bottom prevents instrument to move freely.

Cases are designed to withstand repeated exposure to 134°C steam temperature.



Steam / EO / Plasma Sterilizable

(External)

SK-19-5(2F) 39(L) \* 27(w) \* 4(H)cm

SK-19-5 (2F)

# Sung Kwang Meditech ← €



For Laparoscope.
Silicon bar is placed to prevent laparoscope to move freely.

Cases are designed to withstand repeated exposure to 134°C steam temperature.

Steam / EO / Plasma Sterilizable

(External)

SK-20-1

48(L) \* 9(w) \* 5(H)cm


# **TELESCOPE CASE**



For Telescope.
For Arthroscope
Silicon bar is placed to
prevent laparoscope to move freely.

Cases are designed to withstand repeated exposure to 134°C steam temperature.

Steam / EO / Plasma Sterilizable

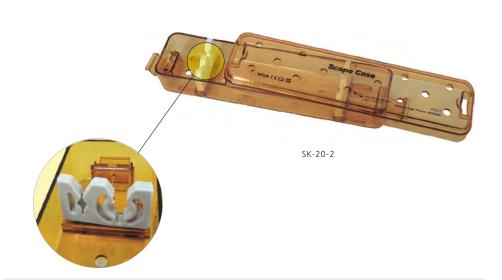


(External)

SK-20-2

28(L) \* 7(w) \* 5(H)cm

# **Operating room**

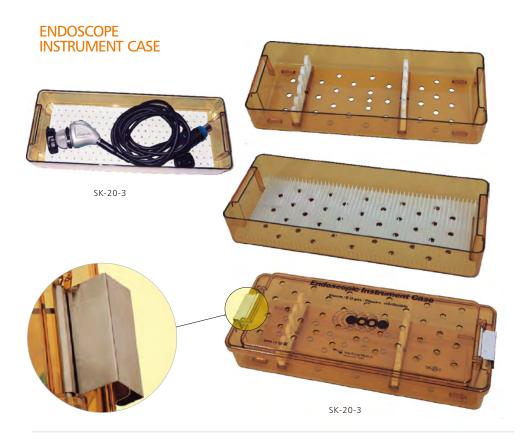

For Endoscope Instrument. Silicon bar is placed to prevent laparoscopic instrument to move freely.

Cases are designed to withstand repeated exposure to 134°C steam temperature.

SK-20-3 is used for endoscope camera case.

Steam / EO / Plasma Sterilizable

(External) 37(L) \* 15(w) \* 7(H) cm


SK-20-3

Silicon bar is placed to prevent laparoscopic instrument to move freely.

Cases are designed to withstand repeated exposure to 134°C steam temperature.

Steam / EO / Plasma Sterilizable

It contains 2pcs or more of laparoscopic instrument.

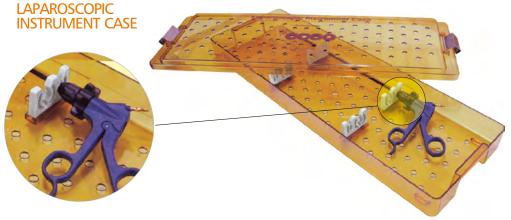



(External)

SK-20-4

60(L) \* 20(w) \* 6(H) cm

# Sung Kwang Meditech ← €



Silicon bar is placed to prevent laparoscopic instrument to move freely.

Cases are designed to withstand repeated exposure to 134°C steam temperature.

Steam / EO / Plasma Sterilizable

It contains 10pcs of laparoscopic instrument.

(External)

SK-21-1

64(L) \* 25(w) \* 15.5(H)cm

Material: stainless steel with marked L/N number, Use: for specimen storage

# **BIOPSY CASE**

R2	Delphain	Pyramida <b>l</b>	L2
R3			L3
R. Jungular	R. Thyroid	L. Thyroid	L. Jungular
R4			L4
R5	R,CN	L,CN	L5
R,Para, T			L,Para, T

SK-21-a-1

Rt level 1		Rt level 1	
Rt level 2		Rt level 2	
Rt level 3		Rt level 3	
Rt deep margin Rt rotter's		Lt deep margin	Lt rotter's
Rt, axillary Rt, sentinal		Lt, axillary	Lt, sentinal

SK-21-a-2



SK-21-a-1(for Thyroid)

1	2	3	5	6	7	4sa	1ab	14
8	9	10	11p	11d	12	14v	450	40

SK-21-a

5	6	8	12a	12b	12c	12h	12p
13	14a	14v	16	16a2	16b1	17	

SK-21-a-4



SK21-a-3 (for stomach)

(External)

SK-21-a for lapa 31.5(W)\*8.5(L)\*3.2(H)
SK-21-a-1 for Thyroid 16.5(W) \* 19(L) \* 3.2(H)

SK-21-a-2 for breast 14.5(W) \* 15(L) \* 3.2(H) SK-21-a-3 for stomach 31.5(W)\*8.5(L)\*3.2(H)

SK-21-a-4 for biliary. Pancreas 30(W)\*8.5(L)\*3.2(H)

# **Operating room**


# POSITION PADS



Special sponge(inner) Good positioner Anti-bedsore effect Very soft, Good elasticity Good endure.

Poly urethane vinyl(outer) Water resistant.


SUPII	NE
Head, f	oot pad
SK-1	22(L) * 16(W) * 11 (5.5)(H)cm
Head p	ad
SK-1-1	22.8(O.D) * 11.5(I.D) * 5(H)cm
Should	er pad
SK-2	23(L) * 14(W) * 7(H)cm

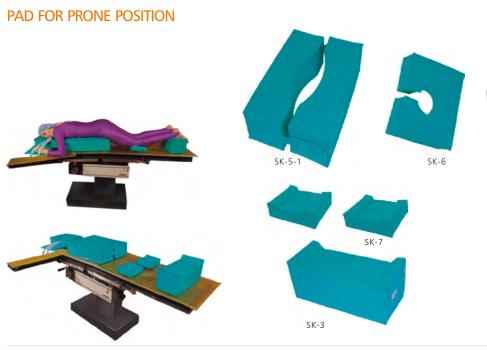



SK-6-1

SK-5-1gel cover

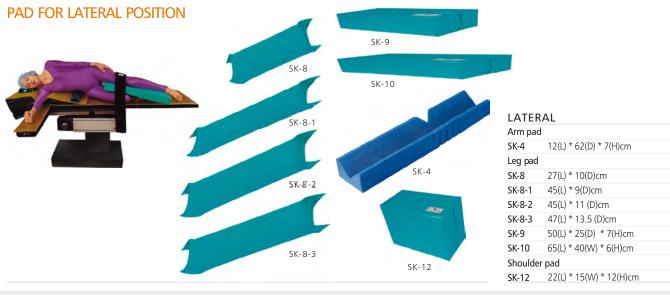

Sung Kwang Meditech Co.,Ltd

www.skmeditech.com

# **GEL PRODUCTS**

Capillary closure is at the root of ulcer formation.

The blockage or damage of capillaries prevents the natural maintenance of healthy skin.

In operation, the use of proper positioning devices prevent excessive pressure build-up and capillary closure.

Evenly distribute weight and provide constant support beneath patients.

Gel products are cost effective, latex free, fast & easy cleaning, radiolucent tissue equivalent.



Sung Kwang Meditech Co.,Ltd

www.skmeditech.com

# **Operating room**



C-DOUGHNUT GEL HEAD PAD



GEL OPHTHALMIC HEADREST



GEL OPHTHALMIC CRADLE HEADREST



GEL HEAD PAD



PRONE GEL HEAD PAD



Doughnut gel head pad provides excellent protection when the patient is positioned for surgery by supine position.

SK-22 (Adult) 20(D) \* 4.5(H) cm SK-22-3 (Child) 15(D) \* 3.5(H) cm SK-22-4 (Baby) 9.2(D) \* 2.4(H) cm

C-doughnut gel head pad provides clear air passage way when the patient is positioned for surgery by lateral or prone position.

SK-22-1 (Adult) 21(D) \* 4.5(H)cm SK-22-1-1 (Child) 15(D) \* 3.5(H)cm SK-22-1-2 (Baby) 9.2(D) \* 2.4(H)cm

It was designed for stable head placement for all ophthalmic and oral or facial procedures.

Also provides patient with comfort while under local anesthesia, minimizing patient movement. It was also designed for surgeon's functional arm support.

SK-22-a (Adult) 31(L)\*27(W)\*10(H) cm SK-22-b (Child) 31(L)\*25(W)\*7.3(H) cm

It was designed for ENT, plastic surgery ,all ophthalmic cases. It helps prevent risk of alopecia especially during long surgery . Cradling of patient's head provides stability .

SK-22-c (Adult) 23(L)\*23(W)\*5.7(H)cm SK-22-d (Child) 19(L) \* 21.6 (W)\*3.8(H) cm

It provides optimum stability of patient's head while in the supine position. Can also be used for position and protect the knee or calf.

SK-22-e (Adult) 18(L)\*18(W)\*8(H) cm SK-22-5 (Baby) 16(L) \* 12 (W)\*1(H) cm

It provides clear air passage way when the patient is positioned for surgery by prone position.

SK-22-2 27(L)\*24(W)\*15(H) cm

# Sung Kwang Meditech C €



Gel body pad provides complete pressure reduction and shares across entire table surface.

SK-23	160(L) * 50(w) * 1.3 (H)cm	
SK-23-1	183(L) * 50(w) * 1.3(H)cm	
SK-23-2	84 (L) * 50(w) * 1.3(H)cm	



This set includes head, torso and foot section pads for use with segment tables.

<b>SK-23-1-1</b> Gel body	pad 3-piece set(a+b+c)
SK-23-1-1a (Head)	50(L) * 23(W) * 1.3(H)cm
SK-23-1-1b (torso)	100(L) * 50(W) * 1.3(H)cm
SK-23-1-1c (foot)	50(L) * 50(W) * 1.3(H)cm

It provides gentle pressure for neo-natal patients in the NIC units.

It is designed for the isolette beds.

SK-23-5	60(L) * 33(W) * 1.3(H)cm
SK-23-5-	1 60(L) * 40(W) * 1.3(H)cm

Gel chair pad is used underneath prominence area in a variety of applications.


SK-24	42(L)*42(W)*1.3(H)cm	
SK-24-2	38(L)*38(W)*1.3(H)cm	
SK-24-2-1	40(L)*40(W)* 1 (H)cm	
SK-24-2-2	40(L)*45(VV)* 1 (H)cm	
SK-24-3	32(L)*32(W)*1.3(H)cm	

Gel embo chair pad is used for iliac crest, under heels and elbows. Is used as kidney rest, axillary roll while patient in lateral position.

Gel embobodypad provides complete pressure reduction and shears cross entiretablesurfa ce.

Gel embo	chair pad
	40(L)*45(W)*1.6 (H)cm
	half body pad
	80(L)*45(W)*1.6 (H)cm
	* * * * * * * * * * * * * * * * * * * *
	body pad
SK-23-3	160(L)*45(W)*1.6 (H)cm



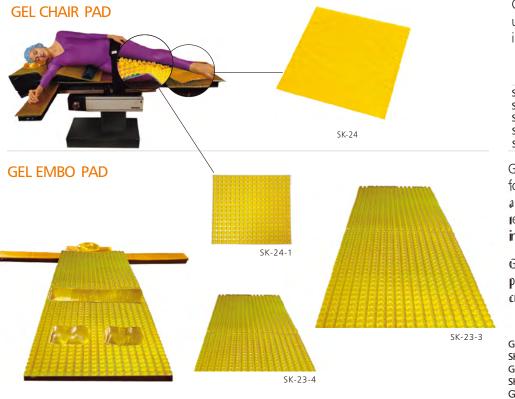

Sung Kwang Meditech Co.,Ltd

www.skmeditech.com

# **Operating room**



Gel knee pad is used for OBGY surgery in lithotomy position. It provides against nerve injury.

SK-25	40(L)*24(W)*1(H)cm
SK-25-1	50(L)*52(W)*1(H)cm

**GEL THYROID PAD** 




Gel lateral pad is designed to accommodate both the shoulder and provides exceptional pressure relief and nerve protection. Gel thyroid pad is designed for thyroidectomy It provides patient with gentle pressure while long surgery with hyperextension position. The tissue friendly feels make it very comfortable. SK-26 70(L) \* 50(W) \* 12(H) / 0.8cm

SK-26-2 45(L)\*30(w)\*7(H) cm

# **GEL PATIENT POSITIONER**

**GEL LATERAL PAD/GEL** 

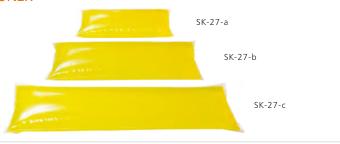

It is good for any procedure. versatile positioner.

SK-27-a (Small)	25(L)*13(W)*3(H)cm
SK-27-b (Large)	35(L)*12(W)*3(H)cm
SK-27-c (Extra long)	46(L)*11(W)*3(H)cm

### **GEL PRONE CHEST PAD**





They were designed for uniformly conform to a patient's shoulder to chest area, with rounded edges for the axilla.

They provide optimum alleviating pressure on the thorax area.

SK-27-e (Left)	70/58(L)*20(W)*5(H)cm
SK-27-d (Right)	70/58(L)*20(W)*5(H)cm

Semi-circular design provides maximum pressure protection while patient is in supine position.

SK-27, 27-4 can be used under the knee, the foot against the board.

While patient is in lateral position, SK-27 can be used for the shoulder protection.

In thyroidectomy, SK-27-4 is used for neck hyperextension.

while patient is in prone position SK-27-1, SK-27-2 can be used under the chest.

SK-27	47(L) * 10(w) * 5 (H)cm	
SK-27-1	47(L) * 15(w) * 8 (H)cm	
SK-27-2	47(L) * 15(w) * 12(H)cm	
SK-27-3	25(L) * 15(w) * 8(H)cm	
SK-27-4	52(L) * 10(w) * 8(H)cm	

# **GEL SHOULDER PAD (CHEST)**



Sung Kwang Meditech Co.,Ltd

www.skmeditech.com

# Sung Kwang Meditech C €



Gel heel pad provides safety and protects heel area body prominence.

SK-28 20(L)\*11(W)\*7(H)cm

It protects heel area without hook and loop attachment. It protects the patient's Achilles tendon area.

It supports heels. Appropriate for recovery .

SK-28-1 16(L)\*8.5(W)\*5(H)cm SK-28-2 17(L)\*18(W)\*1.3(H)cm

It supports heels appropriately for recovery .

SK-28-3 17(L)\*23(W)\*1.3(H)cm

It is applied for supine position for protection of bony prominence area.

SK-29 55(L)\*15(W)\*1.3(H)cm

It is applied for lateral position. Designed to provide arm support and stability.

SK-29-1 (Small) 41(L)\*15(W)\*5(H)cm SK-29-1a (Large) 55(L)\*15(W)\*6(H)cm

It is used for nerve protection in fastening wrist and ankle.

SK-17-a Gel Wrist protector

26(L)\*6(W)\*1(H)cm

SK-17-b Gel Ankle protector 33(L)\*6(W)\*1(H)cm

SK-17-c Gel Armboard binder 21(L)\*6(W)\*1(H)cm

It provides comfort and distributes pressure in prone position spine surgery.

SK-5-1a 52(L) \* 31(W) \* 1.3(H)cm SK-17-1 75(L) \* 24.5(W) \* 1.3(H)cm

It provides comfort and safety in prone position spine surgery.

**SK-17-1-1** 1.3(thick)cm

# **Operating room**

# **GEL FOREARM PROTECTOR**



SK-17-2

**GEL BODY BELT** 



SK-17-3



**GEL SHOULDER PROTECTOR** 



SK-17-3-1

It is used for nerve protection of arm and elbow in supine position surgery.

SK-17-2 Gel forearm protector (L) 46(L)\*21(W)\*0.8(H)cm SK-17-2-1 Gel forearm protector (S)

38(L)\*11.5(W)\*0.8(H)cm

It helps safety, stability and tissue protection during surgery.

SK-17-3 50(L)\*8.5(W)\*0.8(H)cm

It is designed for patient in trendelenburg down position to protect shoulder areas.

16(L)\*15(W)\*2(H)cm

68(L)\*38(W)\*1(H) cm

**GEL ALLEN PAD** 



for nerve protection.

for nerve protection.

SK-17-3-1

SK-17-4-1

**GEL ANKLE HOLDER** (2 PIECE SET)



It is used in lithotomy position

It is used lithotomy position

SK-17-5 40(L)\*3(W)\*1(H)cm

**GEL SNOW PAC** A-1 type COVER



SK-17-6

Designed for A-1 type snow pac. It provides comfort.

**GEL SNOW PAC** B type COVER



100(L)\*73(W)\*1(H)cm SK-17-6

Designed for B type snow pac.

SK-17-6-1

SK-17-6-1 88(L)\*72(W)\*1(H)cm

Designed for C type snow pac.

**GEL SNOW PAC** C type COVER




SK-17-6-2


SK-17-6-2 30(L)\*40(W)\*1(H)cm

Sung Kwang Meditech Co.,Ltd

www.skmeditech.com

# SILICON

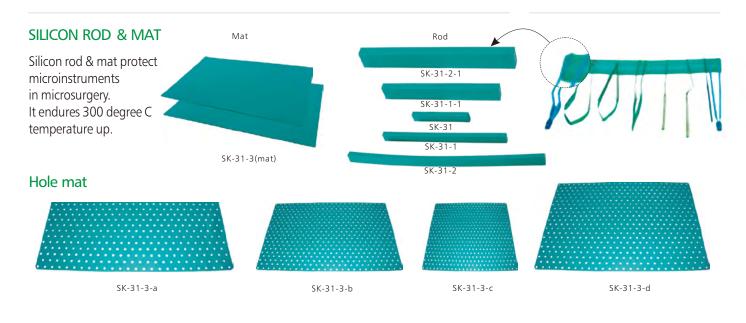

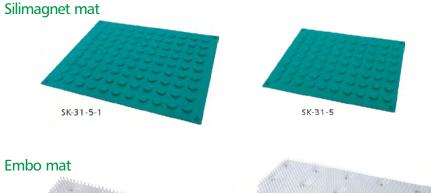



SILICON ROD	
SK-31(Rod)	1.5(L) *10(W) * 0.8(H)cm
SK-31-1(Rod)	1.5(L) * 20(w) * 0.8(H)cm
SK-31-1-1(Rod)	1.5(L) * 20(w) * 2.8(H)cm
SK-31-2(Rod)	1.5(L) * 30(w) * 0.8(H)cm
SK-31-2-1(Rod)	1.5(L) * 30(w) * 2.8(H)cm
SILICON MAT	
SK-31-3(Mat)	50(L) * 33(w) * 0.2(H)cm
SILICON HOLE MAT	Г
SK-31-3-a	50(L) * 25(w) * 0.2(H) cm
SK-31-3-b	33(L) * 25(w) * 0.2(H) cm
SK-31-3-c	25(L) * 25(w) * 0.2(H) cm
SK-31-3-d	50(L) * 33(w) * 0.2(H)cm
SILICON EMBO MA	AT
SK-31-4	23.5(L) * 15 (w) * 1.5(H)cm
SK-31-4-1	37(L) * 24.3 (w) * 1.5(H)cm
SILIMAGNET MAT	
SK-31-5-1	40(L) * 35(w) * 0.6(H)cm
SK-31-5	38(L) * 30(w) * 0.6(H)cm

### **SILICON TUBE**



# SILI SUCTION TUBE (SK-31-4-4)

In neurosurgery , used for suction tube. It provides elasticity, flexibility.

SK-31-4-3 5.0(I.D)\*7.0(O.D) mm SK-31-4-4 6.5(I.D)\*11.0(O.D) mm

# **STERICART**



# 제품특징

- 하중 : 300kg 이상

- 바퀴 : 내열성 ( 200℃), 내부식성

외부 : SUS 304하부 : 충격보호대장치

- 용도 : 병원공급실에서 소독기(Autoclave)

속에 소독품을 실어 카트째 소독할

수 있는 카트



cat.no	SK-42
규 격	93(가로) X 58(세로) X 146(높이)mm

size:(cm)

가로

90

120

150

# **WIRECART**



# 제품특징

- 하중 : 300kg 이상

- 구조 : 무게 견딜 수 있도록 선반이중보강 , 조립식

- 자재 : 철+크롬도금- 바퀴 : 12.5 cm(지름)



Cat.no	폭	가로	높이	Cat.no
SK-40-1	46	90	80	SK-40-4
SK-40-1-1			120	SK-40-4-1
SK-40-1-2			150	SK-40-4-2
SK-40-1-3			180	SK-40-4-3
SK-40-2		120	80	SK-40-5
SK-40-2-1			120	SK-40-5-1
SK-40-2-2			150	SK-40-5-2
SK-40-2-3			180	SK-40-5-3
SK-40-3		150	80	SK-40-6
SK-40-3-1			120	SK-40-6-1
SK-40-3-2			150	SK-40-6-2
SK-40-3-3			180	SK-40-6-3

Sung Kwang Meditech Co.,Ltd

www.skmeditech.com

# SNOW PAC



It is used for lugeric patient and unpositioned handicapped person. Make position easy and fast. Then, suction equipment is used.

RH B type 83(L) \* 108(W)cm HB-102a 33cm

# SEVERAL DEVICES



# Rehabilitation

# GEL WHEELCHAIR CUSHION PAD WITH COVER

It was made with a high-quality polyurethane gel and PU foam.

It provides excellent pressure relief.

The gel material will not migrate, leak, evaporate, harden or dry out if the cushion is accidentally punctured.

The material is highly resistant to bacteria and bodily fluids, and can be easily cleaned with warm soapy water or a disinfectant.


It offers excellent pressure reduction and pain relief combination for people with low to moderate risk of skin breakdown.

It provides excellent immersion stability and pressure distribution.

It provides great comfort and excellent positioning stability. The gel material also helps dissipate heat buildup.


4" cushion provides long-term positioning support and excellent comfort. It is applied for individuals with a history of skin breakdown or those who suffer from pressure sores.

# GEL WHEELCHAIR CUSHION 2"

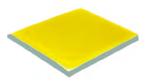

SK-1616-2	16(L)x16(W)x2(H)Inch
SK-1816-2	18(L)x16(W)x2(H)Inch
SK-1818-2	18(I)x18(\M/)x2(H)Inch

# GEL WHEELCHAIR CUSHION 3"

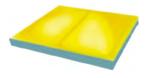


SK-1616-3	16(L)x16(W)x3(H)Inch
SK-1816-3	18(L)x16(W)x3(H)Inch
CV_1Q1Q_2	19/1/v19/\///v3/H\Inch

# GEL WHEELCHAIR CUSHION CURVED 3.5"



SK-C1616-3.5	16(L)x16(W)x3.5(H)Inch
SK-C1816-3.5	18(L)x16(W)x3.5(H)Inch
SK-C1818-3.5	18(L)x18(W)x3.5(H)Inch

# GEL WHEELCHAIR CUSHION MORE CURVED 4"

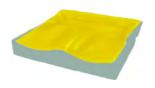

SK-MC1616-4	16(L)x16(W)x4(H)Inch
SK-MC1816-4	18(L)x16(W)x4(H)Inch
SK-MC1818-4	18(L)x18(W)x4(H)Inch

# WHEELCHAIR CUSHION COVER








Bottom cover is made by non slip material.
Side cover is made by net shaped material.
Top cover is made by very soft, skinlike material.

SK-W-cover


# Sung Kwang Meditech C €

<b>GEL WHEELCHAIR</b>	
<b>EMBO CUSHION</b>	2"

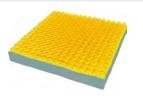


SK-E1616-2	16(L)x16(W)x2(H)Inch
SK-E1816-2	18(L)x16(W)x2(H)Inch
SK-E1818-2	18(L)x18(W)x2(H)Inch

GEL WHEELCHAIR EMBO CUSHION 3"

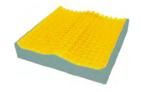


SK-E1616-3	16(L)x16(W)x3(H)Inch
SK-E1816-3	18(L)x16(W)x3(H)Inch
SK-E1818-3	18(L)x18(W)x3(H)Inch

GEL WHEELCHAIR EMBO CURVED CUSHION 3.5"

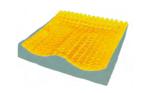

SK-EC1616-3.5	16(L)x16(W)x3.5(H)Inch
SK-EC1816-3.5	18(L)x16(W)x3.5(H)Inch
SK-EC1818-3.5	18(L)x18(W)x3.5(H)Inch

GEL WHEELCHAIR EMBO MORE CURVED CUSHION 4"

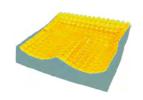

K-EMC1616-4	16(L)x16(W)x4(H)Inch
SK-EMC1816-4	18(L)x16(W)x4(H)Inch
K-EMC1818-4	18(L)x18(W)x4(H)Inch

# GEL PRODUCTS-ACCESSORY

EMBO GEL WHEELCHAIR ARMBOARD PAD





Adhesive tape attached in bottom prevents move freely.

**GEL TOILET COVER** 



It is very effective for users with existing skin breakdown.

23(L)\*7(W)\*1.6(H) cm

RH-4 (Dual)+cover 39(L)\*10(W)\*1(H) cm

**GEL COMMODE SEAT** 





RH-6 (close type) 45(L)\*45(W)cm RH-6-1(open type) 45(L)\*46(W)cm

# Rehabilitation

# GEL WHEELCHAIR BACK CUSHIONS WITH COVER

Polyurethane gel and PU foam helps comfort and durability for long-term.

The embossed gel feature aids comfort for long-lasting back and lumbar support in long-term seating environments.



# GEL HORSE'S EQUESTRIAN PAD

# RECTANGULAR ROUND





SK-EQ-a Rectangular shape with cover 36(L)x50(W)x1(H)cm
SK-EQ-a cover Rectangular shape cover only 36(L)x50(W)x1(H)cm
SK-EQ-b Round shape with cover 40(L)x55(W)x1(H)cm
SK-EQ-b cover Round shape cover only 40(L)x55(W)x1(H)cm

# SILVER PRODUTS

#### GEL EMBO CUSHION PILLOW / MATTRESS

Embo gel cushion pillow with cover (memory foam+embo gel)

### CHARACTERISTICS:

C-spine supporting effect is developed by medical basis. Softly, as humanengineering structure, it supports head. After deep sleeping, it causes comfort.

As urethane foam materials made by physical method.

As urethane foam materials made by physical method, although it is used for a long time, it doesn't transform.



+ embo gel pad Color : silver, ivory Kind : therapeutic, antibacterial, deodorizational.

Outer: 100% cotton velvet with

Materials: slowly recovery

silvercholoid coating

Inner: cotton 75% polyester 25%

polyuréthane foam 100%

RH-1	50(L)*33(W) cm
RH-2 (twin)	RH-2-1 (full)
RH-2-3 (queen)	RH-2-4 (king)
Cover Only	
RH-2-cover	RH-2-1-cover
RH-2-3-cover	RH-2-4-cover

Sung Kwang Meditech Co.,Ltd

www.skmeditech.com

### SEVERAL KINDS OF X-RAY APRON FOR SURGICAL PATIENT





SK-15





SK-15-1





SK-15-4





SK-15-5

Personal protection in the medical use of radiation.

Modern radiology-supported diagnostic and therapeutic procedures extend the wear time for protective clothing.

light weight, high quality, cost effective.

All apron items are available.

inner: natural rubber-based lead rubber flexible wearing comfort and best

outer: nylon fabric, etc satisfies workmanship hygienic requirements and the demands of medical use.

> Pb equivalent: 0.25mm Pb equivalent: 0.35mm Pb equivalent: 0.5mm

#### MALE

	Sm	ıall	M	ed	Lar	ge	X	L
	L	W	L	W	L	W	L	W
Frontal type	90	54	95	56	100	58	105	62
X-band type	91	54	93	56	97	58	101	61
Front/Back type	94/84	51	98/88	52	102/91	54	107/95	56
Vest/Skirt type	57/50	46/40	58/51	48/42	59/52	51/44	60/54	53/46

#### **FEMALE**

	Sm	nall	M	ed	Laı	ge	Х	L
	L	W	L	W	L	W	L	W
Frontal type	82	47	87	50	92	52	97	55
X-band type	82	47	86	50	91	52	95	55
Front/Back type	85/74	43	90/79	46	95/84	49	100/89	52
Vest/Skirt type	53/46	43/37	54/47	45/38	56/49	48/42	57/50	50/44

±1~2Cm, unit : Cm

# **DROP OFF**





Remove your apron quickly and easily with shoulder valcro.

Available to SK-15, SK-15-1, SK-15-4, SK-15-5

# **Lead APRONS**

### SEVERAL KINDS OF X-RAY APRON FOR SURGICAL PATIENT







Thyroid cover protects thyroid gland in the medical use of radiation.

And disposable cover is also helpful in cleanliness.

All apron items are available.







inner: natural rubber-based lead rubber flexible wearing comfort and best

outer: nylon fabric, etc satisfies workmanship hygienic requirements and the demands of medical use.

> Pb equivalent: 0.25mm Pb equivalent: 0.35mm Pb equivalent: 0.5mm









Color type

Dark Blue



SK-602

SK-601

SK-16-3





Like the Page Cap, our Regular Apron Cap is comfortable and lightweight. An adjustable Velcro tab allows custom fitting.

Available in any apron fabric. Please specify when ordering. Order sizes : Male, Female Protection: 0.25mm front

Closures: Velcro

Lead cap(W)	SK-16-3
Lead cap(M)	SK-16-4
Lead gonad shield	SK-16-5 (M,F)

### SEVERAL KINDS OF X-RAY APRON FOR SURGICAL PATIENT





SKN-15





SKN-15-1





SKN-15-4





SKN-15-5

# Personal protection in the medical use of radiation.

Modern radiology-supported diagnostic and therapeutic procedures extend the wear time for protective clothing.

light weight, high quality, cost effective.

All apron items are available.

inner: lead free metal composites flexible wearing comfort and best

outer: nylon fabric, etc satisfies workmanship hygienic requirements and the demands of medical use.

> Pb equivalent: 0.25mm Pb equivalent: 0.35mm Pb equivalent: 0.5mm

### **DROP OFF**





Remove your apron quickly and easily with shoulder valcro.

Available to SKN-15, SKN-15-1, SKN-15-4, SKN-15-5

# Non Lead APRONS

### SEVERAL KINDS OF X-RAY APRON FOR SURGICAL PATIENT













SKN-15-3-2







SKN-16- a

All apron items are available.

inner: lead free metal composites flexible wearing comfort and best

outer: nylon fabric, etc satisfies the workmanship hygienic requirements and the demands of medical use.

> Pb equivalent: 0.25mm Pb equivalent: 0.35mm Pb equivalent: 0.5mm

Thyroid cover protects thyroid gland in the medical use of radiation.

Color type







Dia Chocolate















Melange Purple Melange Orange Melange Blue Melange Navy







Melange Khaki

Black Dot Honeycomb Grey











SKN-16-5 (M)

Like the Page Cap, our Regular Apron Cap is comfortable and lightweight. An adjustable Velcro tab allows custom fitting.

Available in any apron fabric. Please specify when ordering. Order sizes : Male, Female Protection: 0.25mm front

Closures: Velcro

SKN-16-3
SKN-16-4
SKN-16-5 (M, F)



Apron holder will accommodate up to 2-10 aprons.
Unit is ideal if wall space is limited.
It is designed to prevent declining on side in moving.
Solid stainless steel construction prevents rust in long-term use.
High durability under heavy use.
Can be assembled after delivering.

Lead curtain provides separated radiation shielded space by simply drawing.

• protection: 0.5 mmpb

• Order made



Protection for pregnant woman barrier: 0.5 mm pb

Order made

SK-16-6	71(L)*63(W)*140(H)cm
SK-16-6-a	50(L)*63(W)*140(H)cm
SK-16-6-1	25(L)*23(W)cm
SK-16-6-2	25(L)*23(W)cm
SK-16-7-1	57(W)*60(L)/83(H) to 146(H)cm
SK-16-7-3	200(W)x 150(H)cm
SK-16-7-4	100(L)*65(W)*0.25(H)cm

### LEAD BARRIER



Lead-glass viewing panel provides exceptional visibility .

• Protection :1.65 mm pb

• Order made

SK-16-7 lead glass barrier	180(L)x90(W) cm
Window Size:	80(L) x 90(W) cm
SK-16-7-a lead barrier	180(L)x 90(W) cm
SK-16-7-b	180(L)x 90(W) cm
Window Size:	20(L) x 30(W) cm
SK-16-7-c	180(L)x 90(W) cm
Window Size:	180(L) x 90(W) cm

Sung Kwang Meditech Co.,Ltd

www.skmeditech.com

Contact to E-mail: skmeditech@skmeditech.com

Tel: 82-31-855-7548 Fax: 82-31-855-7549

#### 99 ALUMILITE



99AL

Manufactured of lightweight nylon, the 99 Ultralite was designed to be comfortable and provide protection for the entire eye area without the need of an added side shield for lateral exposure.

The curved front and temples with rubber tips assures a close secure fit that also provides some splash protection. An adjustable "snap tip" cord is standard with each model.

- ♦Weight: 70 grams
- ◆Lead Equiv Front: 0.75mm (+/-0.05) Lead glass
- ◆Frame Material: Durable Nylon
- ◆Available In: Plano (non-prescription), single Rx (low power), Bifocal Rx (low power)
- ◆AVAILABLE COLORS: Silver.

ALUMILITE 99AL

A panoramic goggle constructed of impact resistant propionate. The frame front has molded side shielding and is fitted to the brow for increased splash protection. The temple length and frame angle are adjustable for a uniform custom fit.

Model 70S comes with Lead side shields for necessary protection from oblique angle exposure.

- ♦Weight: 80 grams
- ◆Lead Equiv Front: 0.75mm (+/-0.05) Lead glass
- ♦Lead Equiv Side Shields: 0.50mm Lead sheet
- ◆Frame Material: Propionate
- ◆Color: Clear/Blue
- ◆Available In: Plano (non-prescription), Single Rx, Bifocal Rx, Trifocal Rx, and Progressive.

# 70S ASTRO FLEX with SIDE SHIELDS



### **COCOONS FITOVER FLEXFIT**





ASTRO FLEX

70s

Cocoon's Flexfit fitover is a versatile, high quality nylon frame. Featuring Flexfit technology, the fitover's temples can be benched (adjusted) easily for a secure, comfortable fit. The large sized model is ideal for fitting over standard sized prescription eyewear with square or rounded lenses. Blacked out side shields provide a streamlined look as well as 0.50mm lead equivalent protection.

◆CE Marked : Yes

Front Lens: 0.75mm Lead Glass
Side Shields: 0.50mm Pb Strips
Available Prescription: None
Prescription Range: None
Weight: 108g

◆Frame Size : 64-18-150 ◆Color : Black

COCOONS FITOVER

RE-FLEXFIT

# Sung Kwang Meditech ← €

# Lead Masks

### **PANORAMA SHIELD**



400PSM

# PANORAMA MASK



450FPM

### **BREAST SHIELD**



The 400PSM Panorama Shield consists of a Face Formed Radiation Protective Lead Impregnated Acrylic Shield that is suspended from a dual adjustable headpiece. The acrylic shield can be positioned at any angle to the face so as to minimize scatter radiation from passing under the shield by side adjustable dials. The headgear is adjustable in two positions. An adjustable headband that saddles the crown of the head while connecting at a 90  $^{\rm o}$  angle to an adjustable headband that surrounds the head circumference from brow to back. The circumference portion of the headband has an adjustable dial in the rear to increase or decrease the hat size for a comfortable, snug fit. Wipe clean with a soft cloth. Abrasive materials may scratch front lens.

Cannot be sterilized with high temperatures or high pressure

♦Weight: 350 grams

◆Protection: .10mm Lead Acrylic ◆Material: Acrylic/Nylon

PANORAMA SHIELD

400PSN/

The 450FPM Panorama Mask consists of a Face Formed Radiation Protective Lead Impregnated Acrylic Shield that is suspended from an adjustable headpiece. The stationary acrylic shield extends from a headpiece brow connector for a length of 8 inches from top to bottom in front of the entire face. A support strap saddles the crown of the head while connecting at a 90  $^{\circ}$  angle to an adjustable headband that surrounds the head circumference from brow to back. The circumference portion of the headband has an adjustable strap in the rear to increase or decrease the hat size for a comfortable, snug fit. Wipe clean with a soft cloth. Abrasive materials may scratch front lens. Cannot be sterilized with high temperatures or high pressure.

♦Weight: 450 grams

◆Protection: .10mm Lead Acrylic ◆Material: Acrylic/Nylon

PANORAMA MASK

450FPM

For use in spinal imaging. Shield features Velcro panels and adjustable belt for positioning. Neck is lead free for complete spine exposure. Great for scoliosis studies.

Protection: .50mm Closures: Buckle

BREAST SHIELD

APSTOLE

Sung Kwang Meditech Co.,Ltd

www.skmeditech.com

### Lead Free

Our lightest 0.5mm lead-free material

- ◆ Environmentally Safe
- ◆ Equivalent to 0.50mm Pb at 100kVp

### **VEST-SKIRT COMBINATION**





This Vest-Skirt combination is the preferred choice for complete protection and weight distribution. The two piece design shifts 70% of the weight directly to the hips, reducing back fatigue. Skirt features no-slip, quick release belt for added support. Double padded shoulders for maximum comfort. Chest pocket standard. (Pictured with Thyroid Tether.)

- ◆Protection: .50mm front /.25mm back
- ◆Closures: Velcro and Quick- Release Buckle

VEST-SKIRT COMBINATION

AP\/S

# WRAP AROUND BACK RELIEF





Offering full protection, this Back Relief design transfers the weight from the shoulders to the hips with the aid of a 6" bonded elastic belt. (See features of our Back Relief belt.)

This apron features Velcro and buckle closures, as well as a Velcro Breakaway side seam for freedom of movement.

Double padded shoulders for added comfort and a chest pocket are standard features.

- ◆Protection: .50mm front/.25mm back
- ◆Closures: Velcro and Quick -Release Buckle

WRAP AROUND BACK RELIEF

APWABR

Designed to transfer the weight from the shoulders to the hips.

The apron features a unique padded belt system constructed with a bonded elasticized panel.

The stretch proof panel will not lose its shape while providing increased comfort and back support through the life of the apron. Extra thick shoulder pads, chest pocket and shoulder strap are standard. Optional tether on shoulder strap available upon request.

- ◆Protection: .50mm front
- ◆Closures: Velcro and Quick Release Buckle

BACK RELIEF

APBR



**BACK RELIEF** 



Sung Kwang Meditech Co.,Ltd

www.skmeditech.com

# 

### **FLEXBACK**





### VEST-SKIRT FLEXBACK





**BUCKLE APRON** 





FLEXBACK w/SURGICAL DROP SHOULDER





Our Flexback apron is one of the most comfortable models you'll ever wear! Criss-crossing back panels lift the weight off the shoulders and support the lower back. Optional shoulder strap available upon request.

Extra thick shoulder pads are standard. Chest pocket standard. Shown with Thyroid Collar.

◆Protection: .50mm front ◆Closures: Velcro

FLEXBACK APFB

This design combines the benefits of a 2-piece apron with the support of our elasticized flexback apron.

You control the amount of back support by adjusting the inner elasticized quick release belt. Chest pocket standard. Shown with Page Cap.

◆Protection: .50mm front/.25mm back ◆Closures: Velcro and Quick - Release Buckle

VEST-SKIRT FLEXBACK

APVSFR

Quick Release Buckle Apron is user friendly with large support straps criss-crossing the back allowing proper placement.

Belt is adjustable in order to accommodate varying sizes. Chest pocket and shoulder pads are standard.

Optional shoulder strap (as shown on Back Relief APBR above) available upon request.

◆Protection: .50mm front ◆Closures: Quick Release Buckle

**BUCKLE APRON** 

APBA

Two great aprons in one! This is the perfect choice for those who want the convenience of a Surgical Drop Apron and the comfort of our Flexback Apron.

Crisscrossing back panels transfer the weight from the shoulders and distribute it to the hips. Lift the Velcro tabs for quick removal! Extra thick shoulder pads and chest pocket are standard. Shown with Apron Cap.

◆Protection: .50mm front ◆Closures: Velcro

FLEXBACK w/SURGICAL DROP SHOULDER

**APFBSD** 

# **Import**

### **SURGICAL DROP-OFF**





Fluoro done Remove your apron quickly and easily without breaking the sterile field. Extra thick shoulder pads for added comfort. Chest pocket and hanging straps are standard.

- ◆Protection: .50mm front
- ◆Closures: Velcro

SURGICAL DROP-OFF

APSD

# **TIE APRON**





Our most basic apron. Easy to put on. Ties are long enough to fit a large spectrum of sizes. Chest pocket standard.

- ◆Protection: .50mm front
- ◆Closures: Ties

TIE APRON

APEA

Crisscrossing back straps keep the shoulders from slipping off.

Velcro tabs make adjustability quick and easy. Extra thick shoulder pads, for comfort and chest pocket are standard.

# **VELCRO ADJUSTABLE**





- ◆Closures: Velcro

VELCRO ADJUSTABLE

◆Protection: .50mm front

APVA

This apron features an adjustable Quick Release Buckle to secure the back panel and easy ties for the front panel allowing complete adjustability during pregnancy. Velcro release shoulders makes removal easy. Apron also includes an additional .50mm equivalent panel. You can select either the Velcro attachment to the marsupial inner pouch or Velcro belt.

- - ◆Closures: Buckle and Ties

MATERNITY APRON

**APMA** 

# MATERNITY APRON





Sung Kwang Meditech Co.,Ltd

www.skmeditech.com

Contact to E-mail: skmeditech@skmeditech.com Tel: 82-31-855-7548 Fax: 82-31-855-7549

◆Protection: .50mm front / 1.0mm (with added panel) and .25 back

# Lead Gloves and Sleeves

# Sung Kwang Meditech C €

# XGUARD GLOVE - RR2, RR3



FREEGUARD MODEL - 1,2 LEAD FREE, LATEX FREE, POWDER FREE



Our Soft-Touch finish creates a smooth interior reducing skin irritation. Anatomically correct curved fingers reduce hand fatigue. Sizing range from 6 to 9 in half sizes for an unequaled comfortable fit. Attenuates up to 50% more direct beam radiation than other leading products. Thinner fingertips provide enhanced tactile sensitivity.

### **ATTENUATION**

RR-2: 0.012" / 0.30 mm					
Kvp	60	90	120		
Cuff	56%	44%	37%		
Palm	57%	45%	38%		
Finger	57%	56%	38%		

	RR-3: 0.024" / 0.60 mm					
	Kvp	60	90	120		
	Cuff 76% Palm 75%		63%	55%		
			62%	54%		
	Finger	75%	62%	54%		

#### **ATTENUATION**

MODEL 1: 0.009" / 0.22 mm				
Kvp	60	90	120	
Cuff	Cuff 51%		34%	
Palm	55%	44%	37%	
Finger	55%	44%	37%	

MODEL 2 : 0.012" / 0.30 mm					
Kvp	60	90	120		
Cuff	Cuff 62%		43%		
Palm	63%	51%	43%		
Finger	61%	49%	42%		

# **PROGUARD SLEEVE**



**APRON SLEEVE** 



PROGUARD Radiation Reducing Sleeve,

when used in conjunction with the

**PROGUARD Radiation Reducing Procedure Glove** creates an attenuative barrier from your fingertips to your biceps from scatter radiation.

- Inner Nitrile coating reducing proteins, antigens and irritation.
- Lightweight material reduces arm fatigue.
- Extreme flexibility does not restrict movement or dexterity.
- Available in three sizes for a comfortable fit.

ATTENUATION					
PRO <i>GUARD</i>	60 kVp	80 kVp	100 kVp	120 kVp	THICKNESS
Sleeve	55%	45%	35%	31%	0.012" - 0.30mm

PROGUARD SLEEVE SL-S/SL-M/SL-L APRON SLEEVE APSLV

# This optional sleeve is the perfect solution for protecting the upper arm and shoulder. Right or left sides available. Velcro tabs allow easy attachment or removal. Permanent attachment is also available. Features inner Velcro panels for custom girth adjustments. Sleeve can also be retrofit to your existing apron.

#### **BREAST SHIELD**



SURGICAL DROP SHOULDER



THYROID COLLARS



MODEL:[AP-TCO]

MODEL:[AP-TCS]

APSTOLE
APDRSH
APTCO / APTCS

Sung Kwang Meditech Co.,Ltd

www.skmeditech.com

Contact to E-mail: skmeditech@skmeditech.com Tel: 82-31-855-7548 Fax: 82-31-855-7549

33

# Non Lead APRONS Color Swatches

#### NYLON

Unless specified otherwise, nylons come standard 210 denier.



### RIPSTOP




### WIPE-AWAY VINYLS\*




# **ANTI-MICROBIAL POLYS\***


Featuring stain resistant and anti-microbial technology only available from Protech Medical. Our poly material will resist a wide variety of bacteria for several years.

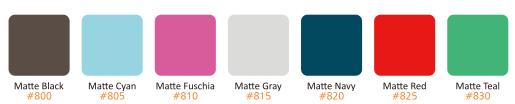


# **PRINTS** | FABRIC COLORS

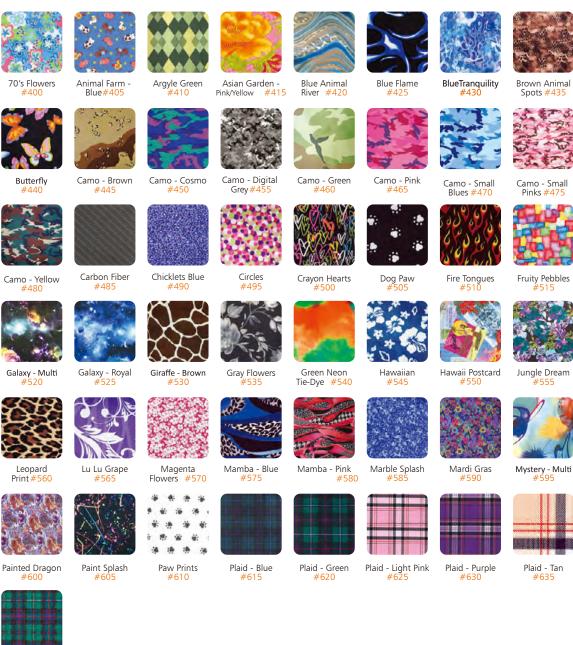

# Plaid - Teal

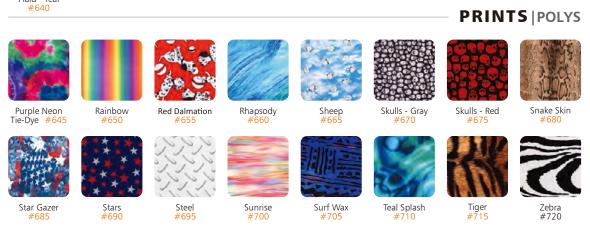

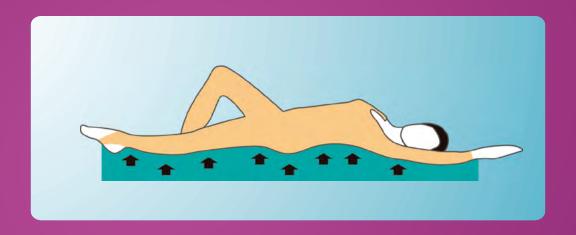



# **Sung Kwang Meditech Co.,Ltd**

Information or ordering call

T) 82-31-855-7548 M) 82-10-4590-7548

F) 82-31-855-7549

É) skmeditech@skmeditech.com

Address) 1014-20, Bukhansan-ro, Jangheung-myeon, Yangju-si, Gyeonggi-do, Korea Sung Kwang Meditech co.,ltd

www.skmeditech.com